



Potency Results

Sample Name: Peach Milano

Client: Azure Remedies

Client Batch ID: Rosin

Pinnacle-Analytics.com

3549 Lear Way, Suite 101

Medford OR 97504

P:(541)300-8217

Sample ID: rC-HS-302-F88

Matrix: Concentrate

Prep Analyst: Jeff A.

Analysis Method: 0668534+1 H4 5-24-2024 #1.lcm

Sampling Method: N/A

Reference Method: JCB 2009: HPLC/DAD

Analysis Batch: 1-8-2025 H4 128, 302, 424, 483 Solids

Date Sampled: 1/8/2025

Date Reported: 1/21/2025


Client License: N/A

3405 172nd St NE Ste 5 PMB #391

Arlington WA 98223

For R&D Purposes Only

Total THC (THCA*0.877+d9-THC) **55.7%**
Total CBD (CBDA*0.877+CBD) **0.34%**
Moisture Content **N/A**

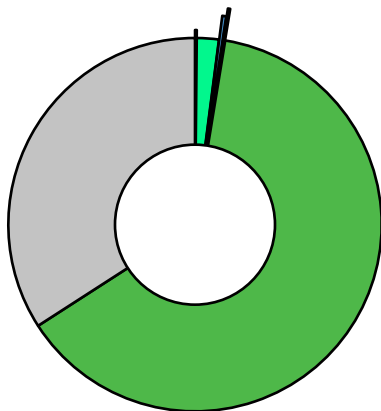

Cannabinoid	% Weight	mg/g
CBDVA	<LOQ	<LOQ
CBDV	<LOQ	<LOQ
CBDA*	<LOQ	<LOQ
CBGA	0.149	1.49
CBG	1.9	19.0
CBD*	0.34	3.4
THCV	<LOQ	<LOQ
CBN	<LOQ	<LOQ
d9-THC*	0.201	2.01
d8-THC*	<LOQ	<LOQ
CBC	<LOQ	<LOQ
THCA*	63.3	633.0
Total Cannabinoids	65.89	659.0

*ORELAP Accredited Analyte

Limit Of Quantitation: 0.09%, analyte not measured

CBGA	d9-THC*
CBG	THCA*
CBD*	Other



These test results may not be altered or reproduced except in full without the permission of Pinnacle Analytics. These results were generated following the Oregon Administrative Rules and in accordance with the NELAP Institute under ORELAP License #4152
 Report generated by Routine_Potency_Rev15_12-23-2024

Kris Ford, PhD
 Lab Director



Quality Control Results

Analyst: Jeff A.

Analysis Batch: 1-8-2025 H4 128, 302, 424, 483 Solids

Pinnacle-Analytics.com
3549 Lear Way, Suite 101
Medford OR 97504
P:(541)300-8217

	Duplicate RPD		LCS % Recovery		Method Blank	
	HS-0-F26-b	Limit	C-SL-010824	Limits	C-SB-010824	Limit
CBDA	3.55%	10%	99.1%	90-110%	<LOQ/2	LOQ/2
CBD	0.546%	10%	99.1%	90-110%	<LOQ/2	LOQ/2
d9-THC	2.61%	10%	103.0%	90-110%	<LOQ/2	LOQ/2
d8-THC	<LOQ%	30%	93.2%	90-110%	<LOQ/2	LOQ/2
THCA	7.96%	10%	98.3%	90-110%	<LOQ/2	LOQ/2

RPD: Relative Percent Difference between unknown sample and its duplicate

LCS: Laboratory Control Sample with known concentration

Case Comments: There were no divergences from ordinary Quality Control procedures or SOPs.



These test results may not be altered or reproduced except in full without the permission of Pinnacle Analytics. These results were generated following the Oregon Administrative Rules and in accordance with the NELAP Institute under ORELAP License #4152
Report generated by Routine_Potency_Rev15_12-23-2024


Kris Ford, PhD
Lab Director